RESMED

 DirectHome
MEDICAL

Presented With Compliments From
www.DirectHomeMedical.com

Swift™ LT

NASAL PILLOWS SYSTEM

User Guide

English • Français • Español • Português

Swift™ LT

NASAL PILLOWS SYSTEM

Rx Only

Thank you for choosing the Swift LT.

Intended Use

The Swift LT channels airflow noninvasively to a patient from a positive airway pressure device such as a continuous positive airway pressure (CPAP) or bilevel system.

The Swift LT is:

- to be used by adult patients (> 66 lb/30 kg) for whom positive airway pressure has been prescribed
- intended for single-patient re-use in the home environment and multipatient re-use in the hospital/institutional environment.

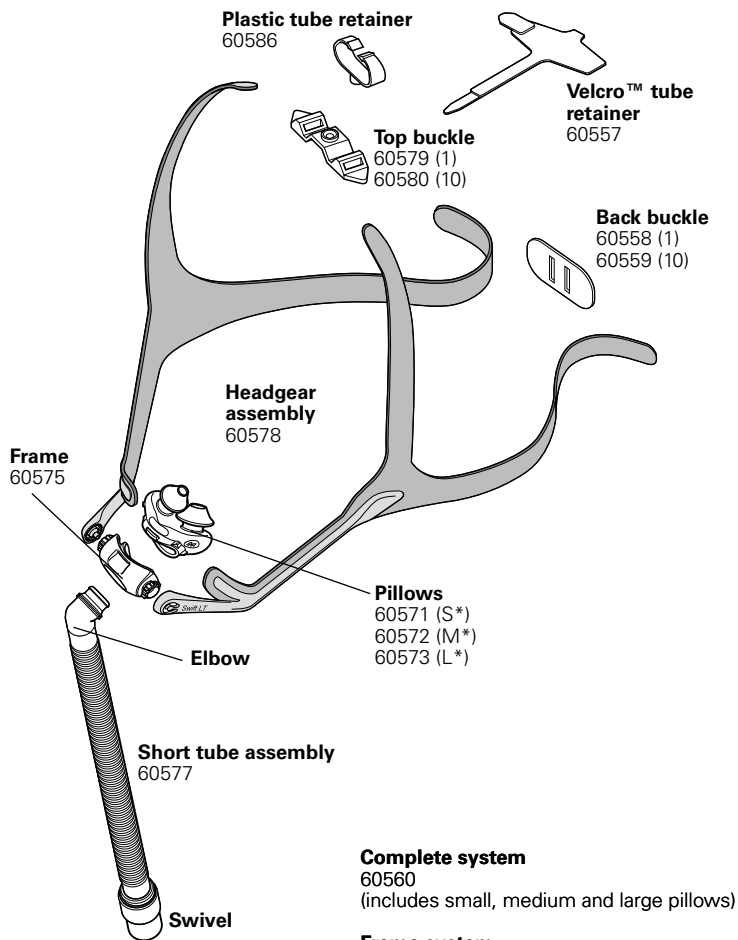
Using the Swift LT

When using the Swift LT with ResMed flow generators that have mask setting options, if available select 'SWIFT'; otherwise select 'MIRAGE' as the mask option.

Notes:

- *The Swift LT is not compatible for use with ResMed AutoSet CS™ 2 and VPAP™ Adapt SV devices.*
- *SmartStop may not operate effectively when using Swift LT with some flow generators that have this feature.*

Parts of the Swift LT & Ordering Information

Complete system

60560

(includes small, medium and large pillows)

Frame system

(pillows, frame, short tube assembly)

S* 60582; M* 60583; L* 60584.

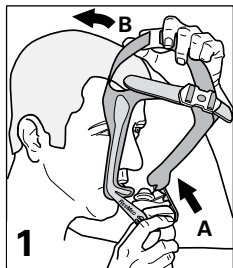
Frame assembly

(frame, short tube assembly)

60581

S* Small; M* Medium; L* Large.


Fitting your Swift LT


(A) Hold the pillows at your nose; **(B)** pull the headgear over your head.

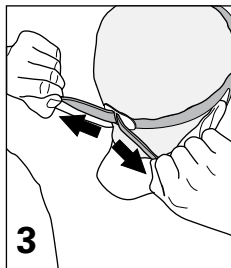
Adjust top straps (overtightening may cause excess pressure on your nose).

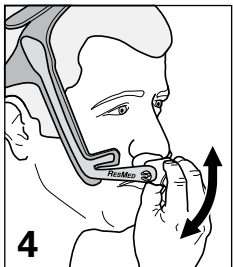Adjust each end of the backstrap (overtightening may cause excess pressure on upper lip).

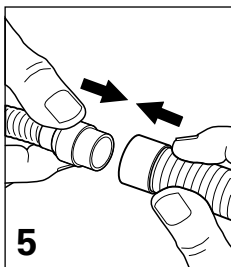Rotate frame to adjust pillow angle to ensure it sits comfortably.

Connect the short tube to the flow generator's air tubing.

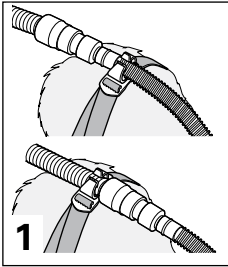
The Swift LT is now fitted and ready for use. Headgear should sit midway between ear and eye.


Using a Tube Retainer (Optional)

You have the choice to use a plastic or Velcro™ tube retainer. Use the version that helps you achieve the most comfortable fit.


Using the Plastic Tube Retainer

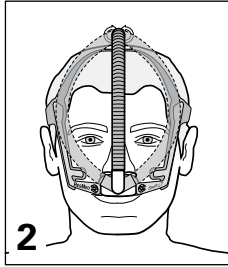Push tube retainer into buckle. Clip either short tube or air tubing to tube retainer.

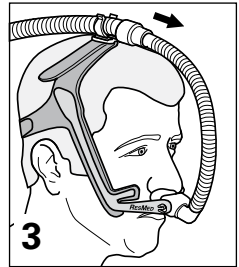
When using the tube retainer, the tubing can be positioned to either side of your face or centered.


If the mask pulls on your nose, increase length of tube between mask and tube retainer.

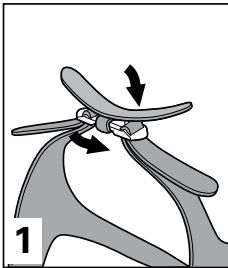
Using the Velcro Tube Retainer

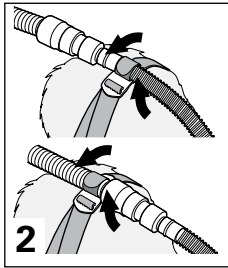Wrap narrow ends of tube retainer around headgear top buckle.

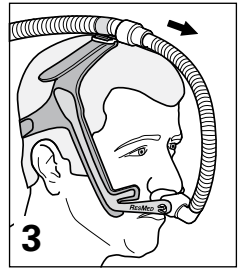
Wrap wide ends of tube retainer around either short tube or air tubing.


If the mask pulls on your nose, increase length of tube between mask and tube retainer.

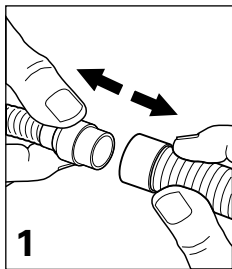
Disassembling your Swift LT


1
Detach air tubing from short tube.

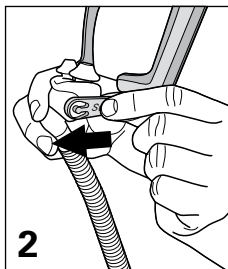2
Remove headgear from frame.

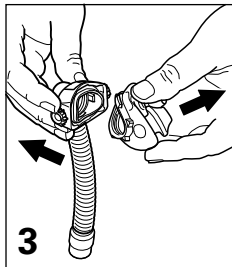3
Remove pillows from frame.

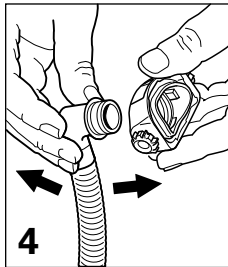4
Remove elbow from frame.

Cleaning your Swift LT in the home

Notes:

- The nasal pillows system and headgear should be cleaned by handwashing.
- The headgear may be washed without being disassembled.

WARNING

Do not use aromatic-based solutions or scented oils (eg, eucalyptus or essential oils), bleach, alcohol or products that smell strongly (eg, citrus) to clean any of the system's components. Residual vapours from these solutions can be inhaled if not rinsed thoroughly. They may also damage the system, causing cracks.

CAUTION

If any visible deterioration of a system component is apparent (cracking, crazing, tears etc.), the component should be discarded and replaced.

Daily/After Each Use

Handwash the separated mask components (excluding headgear) by gently rubbing in warm (approx. 86°F/30°C) water using mild soap for one minute. Rinse all components well with drinking quality water and allow them to air dry out of direct sunlight.

Weekly

Handwash the headgear and tube retainer in warm (approx. 86°F/30°C) water using mild soap. Rinse well and allow them to air dry out of direct sunlight.

Reprocessing your Swift LT between Patients

The Nasal Pillows System should be reprocessed when used between patients. Cleaning, disinfection and sterilization instructions are available from the ResMed website, www.resmed.com/masks/sterilization/americas.

If you do not have Internet access, please contact your ResMed representative.

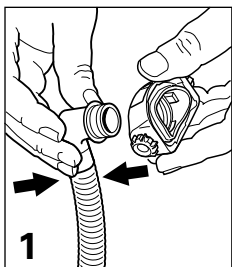
Storage

Ensure that your Swift LT is thoroughly clean and dry before storing it for any length of time. Store the mask in a cool dry place out of direct sunlight.


Disposal

The Swift LT does not contain any hazardous substances and may be disposed of with your normal household refuse.


Reassembling your Swift LT

1 Insert elbow into frame.

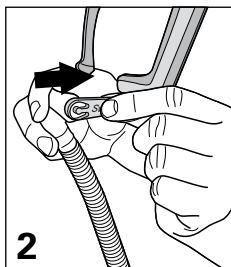
2 Slide headgear clips onto frame.

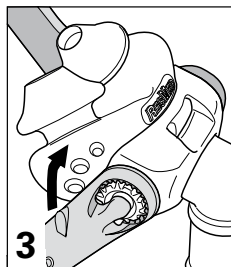
3 Insert pillows into frame, ensuring positioning markers on cushion and headgear are on the same side.

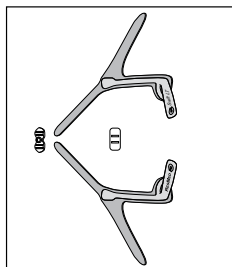
Reassembling the headgear

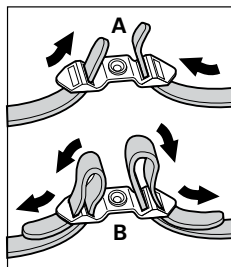Lay out the headgear as shown, writing facing upwards.

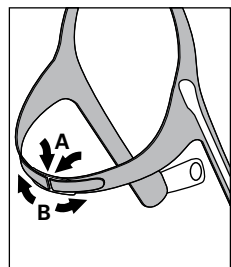(A) Thread the top straps through the top buckle, then (B) back over to secure.

(A) Thread both backstraps through the back buckle (Velcro™ facing outwards), then (B) fold back over to secure.

Limited Warranty

ResMed warrants that your ResMed system, including its components: frame, pillow sleeve, headgear, tubing and other accessories ("Product") shall be free from defects in material and workmanship for a period of ninety (90) days from the date of purchase by the initial consumer.

For further details refer to the Warranty.

Troubleshooting

Pillows won't seal properly, are uncomfortable or cause red marks.

Possible reason	Solution
Pillows may have been fitted incorrectly, incorrectly adjusted, or the headgear is too tight.	Carefully follow instructions in "Fitting your Swift LT". Check headgear is not over-tightened. Check ResMed logo on top of cushions is facing outwards.
Pillows are wrong size.	Talk to your clinician.
Pillows may be dirty.	Clean pillows according to instructions in "Cleaning your Swift LT in the Home".

Nasal pillows system is too noisy

System is assembled incorrectly.	Reassemble the Swift LT. Make sure the cushion is fully inserted into the frame.
----------------------------------	--

Before using your Swift LT

WARNINGS

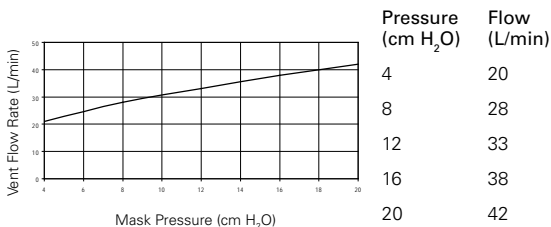
- The vent holes must be kept clear.
- The mask should not be used unless the CPAP system is turned on and operating properly.
- Follow all precautions when using supplemental oxygen.
- Oxygen flow must be turned off when the flow generator is not operating, so that unused oxygen does not accumulate within the flow generator enclosure and create a risk of fire.
- At a fixed flow rate of supplemental oxygen flow, the inhaled oxygen concentration varies, depending on the pressure settings, patient breathing pattern, mask, point of application and leak rate.
- The technical specifications of the mask are provided for your clinician to check that they are compatible with the flow generator. If used outside specification or if used with incompatible devices, the seal and comfort of the mask may not be effective, optimum therapy may not be achieved, and leak, or variation in the rate of leak, may affect the flow generator function.
- Discontinue using the Swift LT if you have ANY adverse reaction to the use of the mask, and consult your physician or sleep therapist.
- As with all masks, some rebreathing may occur at low CPAP pressures.
- Refer to your flow generator manual for details on settings and operational information.
- Remove all packaging before using the mask.

Technical specifications

Note: The mask system does not contain latex, PVC or DEHP materials.

Pressure-Flow Curve

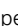
Note: The mask contains passive venting that meets the requirements for protection against rebreathing as specified in ISO 17510-2. As a result of manufacturing variations, the vent flow rate may vary.

Dead Space Information	Physical dead space is the empty volume of the mask to the end of the swivel. Using the large cushions this is 94 mL.
Standards	CE designation in accordance with EC directive 93/42/EEC, class IIa. ISO 14971, ISO 13485, ISO 10993-1 and ISO 17510-2.
Therapy Pressure	4 to 20 cm H ₂ O.
Resistance	Drop in pressure measured (nominal) at 50 L/min: 0.54 cm H ₂ O and at 100 L/min: 1.95 cm H ₂ O.
Sound	DECLARED DUAL-NUMBER NOISE EMISSION VALUES in accordance with ISO 4871. The A-weighted sound power level of the mask is 25 dBA, with uncertainty 3 dBA. The A-weighted sound pressure level of the mask at a distance of 1 m is 17 dBA, with uncertainty 3 dBA.
Environmental Conditions	Operating temperature: +41°F to 104°F (+5°C to +40°C) Operating humidity: 10%–95% non-condensing Storage and transport temperature: -4°F to 140°F (-20°C to +60°C) Storage and transport humidity: 10-95% non-condensing.
Gross Dimensions	380 mm (H) x 57 mm (W) x 62 mm (D). Mask fully assembled with short tube assembly – no headgear.

Note: The manufacturer reserves the right to change these specifications without notice.

Symbols

⚠ Caution, consult accompanying documents; LOT Lot number; REF Part number;  Temperature limitation;  Humidity limitation;  Does not contain latex;  Indicates a Warning or Caution and alerts you to a possible injury or explains special measures for the safe and effective use of the device.

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ResMed Ltd (Manufacturer) 1 Elizabeth Macarthur Drive Bella Vista NSW 2153 Australia
ResMed Corp (US Designated Agent) 14040 Danielson Street Poway CA 92064-6857 USA
ResMed (UK) Ltd (EU Authorized Representative) 96 Milton Park Abingdon Oxfordshire OX14 4RY UK

ResMed Offices Australia, Austria, Brazil, China, Finland, France, Germany, Hong Kong, India, Japan, Malaysia, Netherlands, New Zealand, Norway, Singapore, Spain, Sweden, Switzerland, UK, USA (see www.resmed.com for contact details).

Swift LT Nasal Pillows System

Protected by patents: AU 785376, DE 20321468.4, DE 20321469.2, DE 20321470.6, DE 20321471.4, DE 20321472.2, EP 1314445, NZ 553756, NZ 553757, NZ 553822, NZ 553824, NZ 553825, US 6581594, US 6823865, US 7159587, US 7318437. Other patents pending.

Protected by design registrations: AU 301943, CN 200430102469.3, JP 1243026, NZ 407161, US D552731. Others designs pending.

Swift, AutoSet CS and VPAP are trademarks of ResMed Ltd and are registered in U.S. Patent and Trademark Office.

Velcro is a registered trademark of Velcro Industries B.V.

© 2008 ResMed Ltd.

608226/1 08 03

SWIFT LT

USER

AMER MULTI