



CHERISH. EVERY. BREATH.

Drive DeVilbiss 5-Liter Oxygen Concentrator (525)

Built upon a legacy of reliable, field-proven concentrator models, the DeVilbiss 5 Liter Oxygen Concentrators are designed with enhanced quality, durability, performance and simplicity in mind.

Assembled in the USA, the DeVilbiss 5 liter concentrators are robust, state-of-the-art oxygen systems that **offer the right balance of cool and quiet performance.**

With improved reliability and reduced operational costs, the 5 liter concentrator is perfect for your patients and the right investment for your business.



CELEBRATING
130
YEARS


For more than 130 years, DeVilbiss products have paved the way for the respiratory industry. Here's another product they'll love from the name you can trust.



5 Liter
Concentrator



10 Liter
Concentrator

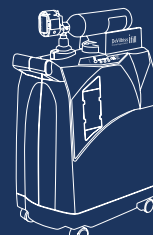

iFill® Personal
Oxygen Station



iGo® Portable
Oxygen
Concentrator



iGo®2 Portable
Oxygen
Concentrator

Performance and Design Features

- Exclusive DeVilbiss OSD® (Oxygen Sensing Device), standard on every unit, ensures patient safety and reliability for longer service intervals
- Readily accessible patient controls, protected cannula fitting and recessed humidifier nook to prevent damage
- Lockable flow meter (Low output flow meter available)
- Convenient top and side handles for easy transport
- Accommodates up to a 50-foot oxygen tubing plus 7-foot cannula (15.3 meters)
- Robust design for improved durability
 - Unit tested in high temperature and high humidity for extended periods of time
 - Proprietary sieve mixture, fill process and assembly
 - Proven rotary valve design
 - Exceptional compressor durability
- Three-year provider service interval
- Auxiliary Oxygen Port



515LF-607 Low Output Flow Meter

Safety Features

- To reduce service and maintenance, alarms do not require a battery. Alarms indicate:
 - Power failure
 - High gas temperature
 - High pressure
 - Low/high flow
 - Low oxygen
 - Service required
- Pressure-compensated flow meter to ensure accurate flow display
- Pressure-relief valve and thermal protection on the compressor
- Flame-retardant cabinet
- Double insulated unit and two-prong plug

Specifications	525DS
Outlet Pressure	8.5 ± 0.5 psig (58.6 ± 3.5 kPa)
Oxygen Percentage	1–5 LPM=87-96%
Power Consumption	310 Watt Average; 275 Watts @ 1.2 LPM & below
Weight	36 lbs. (16.3 Kilograms)
Dimensions	24.5" H x 13.5" W x 12" D (62.2 x 34.2 x 30.4 cm)
Warranty	3 years limited
Alarms	Audible and visual high pressure, low/high flow, low oxygen, high gas temperature, power failure, Oxygen Sensing Device (OSD)

Ordering Information	
525DS	5-liter oxygen concentrator
525DS-Q	Ultra-quiet 5-liter oxygen concentrator



For more than 130 years, DeVilbiss products have paved the way for the respiratory industry. Here's another product they'll love from the name you can trust.

