



Performance
Made Compact



[CLICK TO LEARN MORE](#)

The next generation of Dynamic Support Technology™


The F&P Evora™ Full Mask is uniquely designed to optimize performance throughout the night.

Dynamic Support Technology combines Stability Wings with the Floating Seal to allow for freedom of movement while keeping the mask in place.

Stability Wings
Offer support to the seal

Floating Seal
Soft, pliable silicone seal decoupled from the frame

Stability

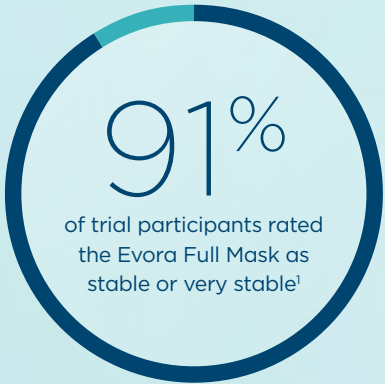

Seal Performance





Secure Comfort

The Floating Seal accommodates for movement while the soft silicone gently wraps under the nose for a comfortable night's sleep.


Comfort






Removable Frame Clip

Enables the headgear settings to stay in place during fitting and cleaning

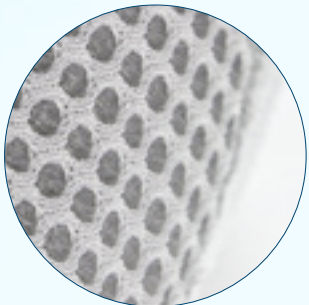

Headgear with VentiCool™ Technology

Breathable headgear designed to allow for heat and moisture to escape

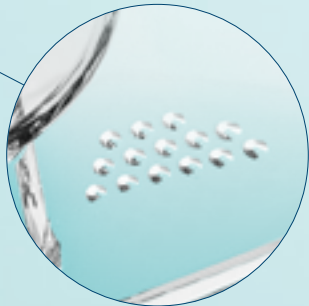

Exhaust Holes

Designed for minimal noise

Sizing Options



XS

Extra Small

S-M

Small to Medium

L

Large



Std.

Standard

XL

Extra Large

Fitting Summary



Guide the headgear
over your head



Hook both the headgear
clips onto the frame



Adjust top, then
bottom straps

For further instructions, please refer to Use and Care Guide.

Product Codes

Code	Description
EVFIXA	Evora Full Mask Extra Small A Model
EVFIMA	Evora Full Mask Small-Medium A Model
EVFILA	Evora Full Mask Large A Model
EVFIXMLA	Evora Full Mask Fitpack A Model

Spare Parts

Universal frame fits all seal sizes

Code	Description
400EVF114	Evora Full Seal - Extra Small Spare
400EVF115	Evora Full Seal - Small to Medium Spare
400EVF116	Evora Full Seal - Large Spare
400EVF121	Evora Full Headgear - Standard Spare
400EVF122	Evora Full Headgear - Extra Large Spare
400EVF131	Evora Full Mask - No Headgear, Extra Small Seal
400EVF132	Evora Full Mask - No Headgear, Small to Medium Seal
400EVF133	Evora Full Mask - No Headgear, Large Seal
400EVF141	Evora Full Headgear Clips & Frame Clip Spare
400EVF151	Evora Full Tube Spare
400EVF161	Evora Full Frame Spare
400FPH341	F&P Swivel White Spare

F&P Evora Full is available
on the F&P myMask™ App

Instructs patients on mask use through step-by-step fitting,
fine-tuning and cleaning videos. Download the F&P myMask App
or scan the QR code to download directly.



1. 40 out of 44 participants rated the Evora Full stability as either very stable or stable. Internal validation study conducted with 44 participants in New Zealand 2020. 2. 39 out of 43 participants rated the overall seal performance of the Evora Full mask as very good or good. Internal validation study conducted with 44 participants in New Zealand 2020. 3. 41 out of 44 rated the Evora Full mask as comfortable or very comfortable. Internal validation study conducted with 44 participants in New Zealand 2020.

Apple, the Apple logo and iPhone are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple Inc., registered in the U.S. and other countries. Google Play and Google Play logo are trademarks of Google LLC. F&P, Evora, myMask, Dynamic Support Technology and VentiCool are trademarks of Fisher & Paykel Healthcare Limited. For patent information, please see www.fphcare.com/ip. REF 625824 REV A 2022-03 © 2022 Fisher & Paykel Healthcare Limited.